ΔZZΔ mumatics™

Certification expertise

Sanitary certification

ASCO Numatics Sanitary diaphragm valves are certified to the European Hygienic Equipment Design Group (EHEDG) document no 8. They are also certified to 3A sanitary standards and approved by the USA Food and drug Administration Agency, FDA.







Functional Safety

551-552-553 pilot valves are certified according to IEC 61508. Functional safety data is available on request. The pilot valves are approved for use up to SIL 3 loop.





ATEX-IECEX certification

Compliance with ATEX directive 94/9/EC for use in potentially explosive atmospheres. Many operators



ASCO Numatics valve islands with G3 electronics support many of the standard communications protocols, including Foundation Fieldbus, DeviceNet, EtherNet/IP, PROFIBUS-DP, PROFINET and MODBUS TCP. There are also many options for ASinterface communications, including solenoid valves, proportional valves and travel feedback.













available with different categories and protection mode (Ex d,Ex ia, Ex m, Ex e mb ...).

Biotechnology& Pharmaceutical









Our expertise, your efficiency

As a company with considerable expertise in providing fluid automation solutions to the process industries, we know that making the right choices can enable efficiency to be improved or maintained.

We offer a broad range of high quality and industry proven fluid control and fluid power solutions and products, including many that are certified to, and meet the needs of the Pharmaceutical and Biotechnology industries

The right fit for your process

With a reputation for quality and reliability our fluid automation solutions are ideal for all critical Pharmaceutical and Biotechnology plant areas. These include the processing of raw materials and finished products as well as the control of utilities.

Processing Raw material delivery to reactors and fermenters, sampling, injection and dosing, filling, coating & packaging - both in dry and liquid forms.

Utilities Pure water generation, CIP-SIP circuit, steam & water sampling, filtration and purification, fluid networks (water, air, nitrogen steam and oxygen), wastewater treatment.



The right solutions for your applications

Our products and solutions are suitable for most Pharmaceutical and Biotechnology applications and we have the expertise to apply them, enabling the efficiency of your plant and processes.

Hygienic Perfectly designed to suite CIP and SIP processes, our aseptic valves are self draining and contain no dead legs or retention areas to reduce bacteriologic risks.

Certification Our Biotechnology and Pharmaceutical products meet mandatory certification requirements for critical, safe and hygienic processes.


Corrosion Resistance Where required, our products are resistant to corrosive and chemical attacks caused by the process and cleaning fluids used in many Biotechnology and Pharmaceutical plants. The materials used will not leach out and adulterate the flowing media.

Diagnostics For easier start-up, commissioning and trouble-shooting, our valves, valve islands and piloting solutions are equipped with sensors, controls and clear diagnostics.

Solutions For ease of installation we offer complete solutions that meet you application needs. We will mount our products in cabinets, on panels, on manifolds or manufacture complete sub systems.





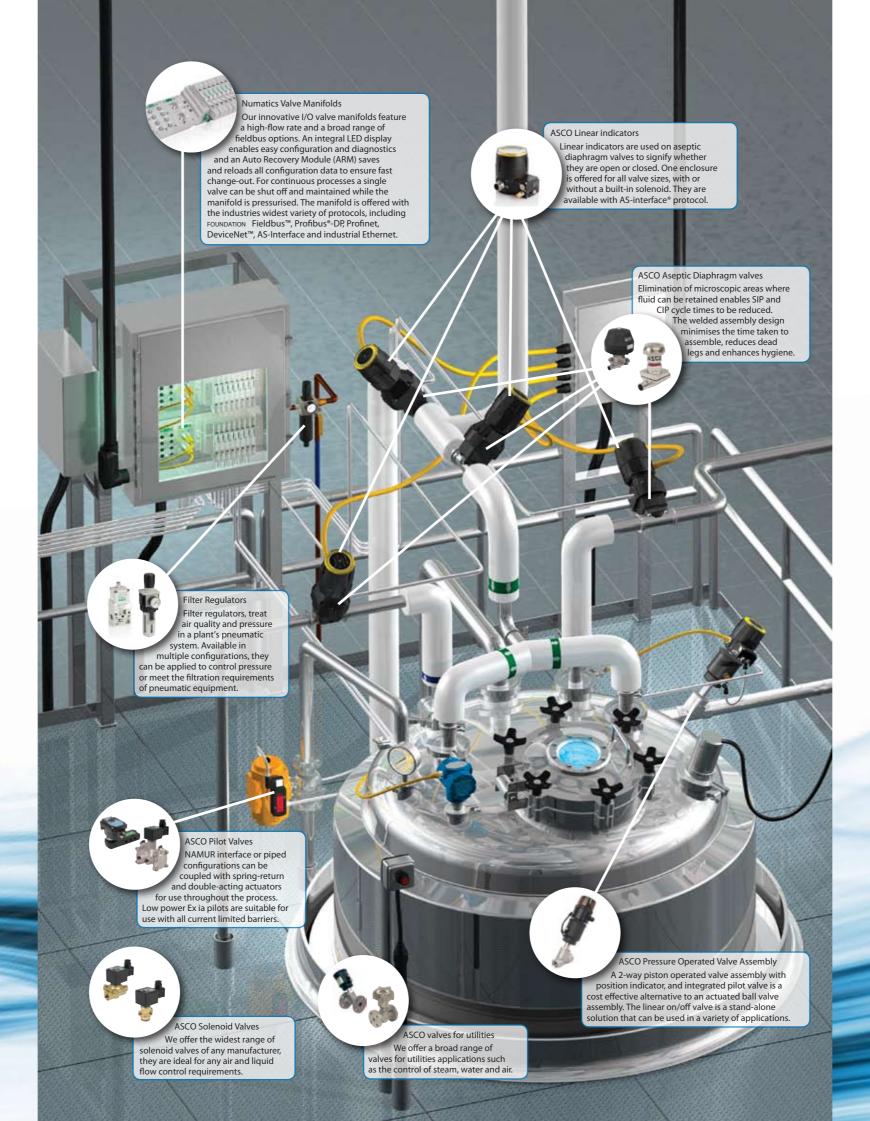



CIP & SIP systems

Our aseptic valves are offered with a range of diaphragm materials and have contoured, polished and self draining bodies.
To ease installation issues, the valves can be provided as a custom



solution with fully assembled and welded inlet and outlet piping. The high quality valves are resistant to steam, acids and other fluids used in typical cleaning or sterilisation systems. In order to meet most common requirements they can be manually operated or operated by air supplied by a solenoid valve or valve island. Our valve islands are available with most industry standard fieldbus protocols.

Valve piloting solutions

Process valves can be controlled by a single pilot valve or by a valve installed on a manifold. Valve manifolds can be assembled directly on the machine or integrated into a specially designed



cabinet that includes all electrical and pneumatic connections. Our valve piloting solutions can be installed in safe areas or in explosive atmospheres.

Pure water generation

This process requires high quality valves, filter elements and flow meters. We have the range, the expertise and the technology to fulfil these functions.



Utilities

Fluids such as steam, dry air, water and gases such as nitrogen oxygen and Carbon Dioxide, can be supplied to the process via ASCO Numatics solenoid or air-operate valves.

